## Introduction

The design specifications and premium quality materials used to make the Shaker Medic Plus ensure the device is capable of being decontaminated using various methods. The Shaker Medic Plus has undergone laboratory testing to prove that being cleaned with enzymatic detergent alone will reduce the biological risk to an acceptable level. If greater sterilisation is required then Shaker Medic Plus can be autoclaved to 134°C / 273 °F at 1Kg/m$^2$ pressure

**For single patient use - See the instruction in the user's manual.**


**For multi-patient use:**

## Method One - Hand Cleansing

1. Disassemble the device by unscrewing the top and tipping out the ball and cone.
2. Fully immerse the parts in an enzymatic cleaner such as Lancerzyme in a clean sink and wash thoroughly using either a soft cloth or cleaning pad to remove any deposits.
3. Rinse thoroughly with deionised or distilled water a minimum of 3 times.
4. Allow to dry naturally. Do not use a heater/hair dryer or other simulated heat source
5. Reassemble and place in a clean container.

## Method Two – Automatic Dishwasher

1. Disassemble the device by unscrewing the top and tipping out the ball and cone.
2. Place all parts in a suitable basket or on the top shelf and wash in an automatic dishwasher. The dishwasher cycle must have the following key elements
   a. Pre-wash - approximately 2 mins
   b. 3 minute detergents wash at 93°C / 199 °F
   c. A complete rinse cycle
   d. Drying stage sufficient to ensure the parts are completely dry.
3. Reassemble and place in a clean container.
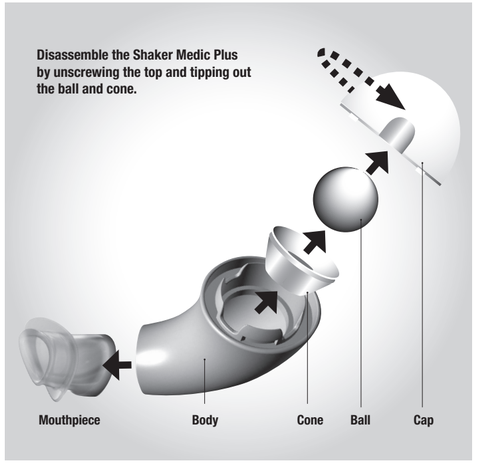
## Method Three - Autoclave


After cleaning method one or method two, transfer the parts to an autoclave and sterilise at 134°C / 273 °F at 1Kg m$^2$ pressure and then place in a sterile container.

**Shaker Medic Plus:** the manufacturer guarantees the decontamination (efficiency of the sterilization) of the product, without any alteration or loss of function, for up to ten (10) processes of autoclaving, except the PVC mouthpiece guaranteed up to three (3) processes of autoclaving. However, the mouthpiece coloration might change after the first autoclaving session. The discolouration itself is not a problem, but make sure that the mouthpiece functionality has not been compromised.

**Note: Shaker Medic Plus cannot be cleaned in a microwave**

## Cleaning and Disinfection

1. Disassemble the Shaker Medic Plus and inspect for signs of damage or wear, if any is evident replace.
2. Prepare a solution of detergent in accordance with the manufacturer instructions in a container large enough to totally submerge the components.
3. Leave the Shaker Medic Plus submerged in the solution for recommend time. Wash thoroughly using either a soft cloth or cleaning pad, ensure all mucus deposits have been removed.
4. Rinse as recommended. Rinse and dry as recommended.
5. Pour a quantity of the chosen disinfectant into suitable container.
6. Leave the Shaker Medic Plus completely submerged in the solution for the recommend time, ensuring that there are no air bubbles.
7. Rinse as recommended - Dry naturally and reassemble.
   **(Do not use hot air or a drying cupboard).**
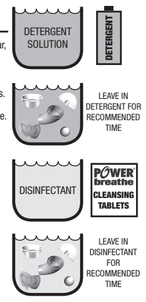
## Disinfectants

To the best of our knowledge the following disinfectants can be used with the Shaker Medic Plus at the stated concentrations only.

| Chemical type | Examples | Solution strength | Comments |
|---|---|---|---|
| Chlorine dioxide Generator | Tristel one day | 20ml in 1 litre of water | Safety data sheet and further information available from www.tristel.com |
| Sodium hypochlorite (NaOCL) | Milton | 1000 ppm | Ensure thorough rinsing, corrosion of the metal parts 4 as will occur if exposed to chlorine for long periods. |
| Sodium dichloroisocyanurate (NaDCC) | POWERbreathe Cleansing Tablets or Presept, Actichlor Sanichlor, Haz-Tab or similar | 1000 ppm | Follow instructions provided by the manufacturer. |

- POWERbreathe International Ltd accepts no liability for damage caused to products if the above procedure and recommended solutions are not used.
- It is the user's responsibility to choose which of the recommended solutions are used within their institution or hospital and we stress that the infection control nurse/department should be consulted when making the choice.
- It is the responsibility of the user to keep them up to date with the latest information from the relevant disinfectant manufacturer concerning instructions, effects, concentrations, and immersion times.
- If your preferred cleaner/disinfectant is not on the recommended list please contact our customer service department Tel:+44 (0) 1926 816100. email enquiries@powerbreathe.com

**References:** Chemical disinfection in hospitals – PHLS, Sterilization, disinfection and cleaning of medical device equipment (MAC Manual) – MHRA

**Note:**

1. This method of cleaning has been independently tested and validated to ensure that if followed the Shaker Medic Plus will be sufficiently clean for re-use.
2. Lancerzyme® is a registered trade mark of Lancer UK Ltd (www.lancer.co.uk) other options may be available locally in other countries
3. The amount of time the device is left in the solution will depend on the disinfectant chosen.
4. The ball in the Shaker Medic Plus is nickel plated carbon steel; in tests some balls have shown signs of discolouration and scratches after prolonged submersion in some cleaning solutions. However, the functionality of the product has not been affected.
5. It is the responsibility of the user to ensure all sterility verifications

**The Shaker Medic Plus cannot be sterilized with ethylene oxide or hot air.**

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This Product is Protected by Patent and Design Copyright ©